


REVISION SCHEDULE		
REVISION #	DESCRIPTION	PLANNER/DATE

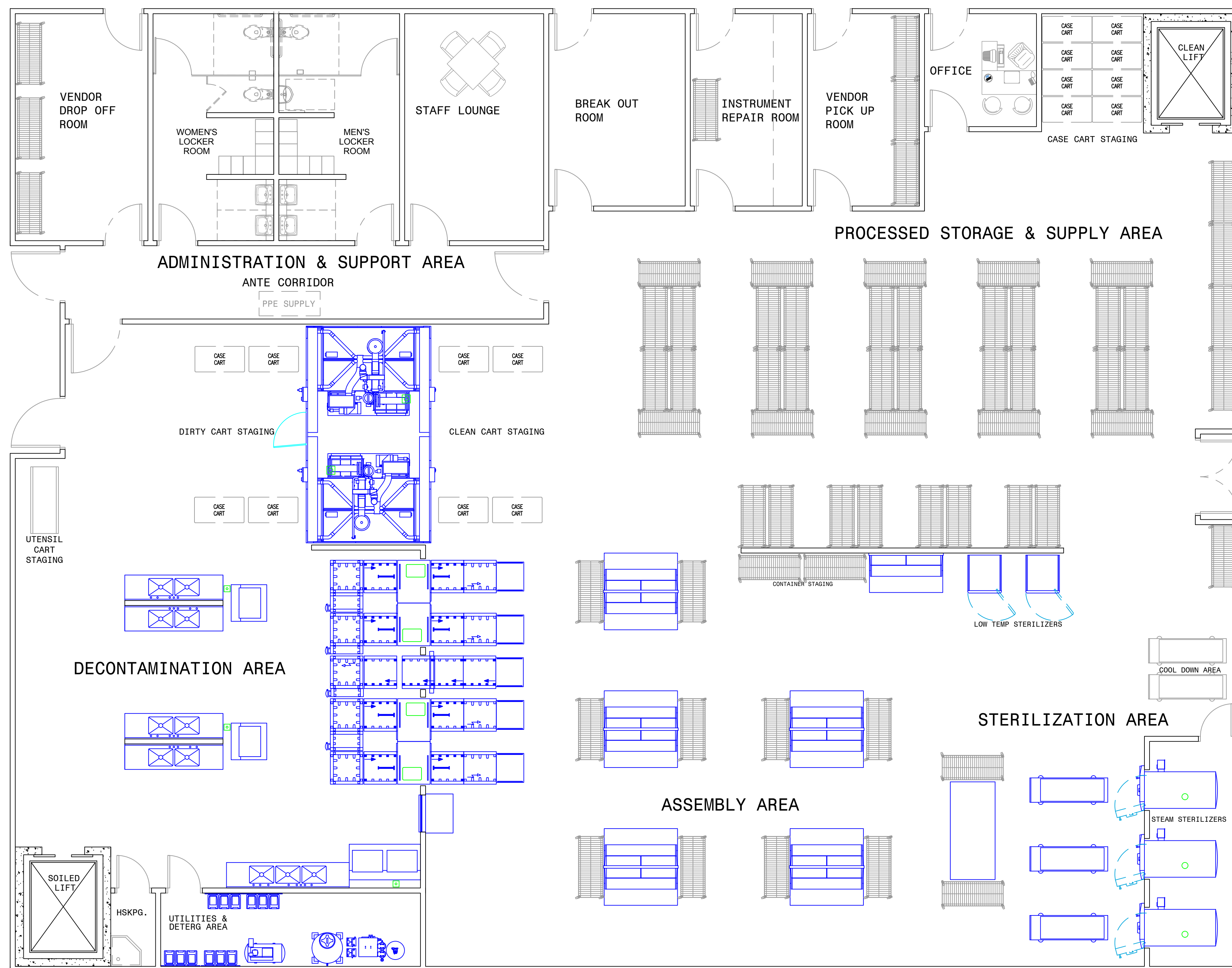


## STERILE PROCESSING DEPARTMENT FLOOR PLAN

**\*THIS SAMPLE LAYOUT IS INTENDED TO ENCOURAGE DISCUSSIONS FOR AN INTEGRATED DESIGN PROCESS. FOR A COMPREHENSIVE CAPACITY ASSESSMENT AND OPTIMIZED DESIGN FOR YOUR PROCESSING DEPARTMENT, PLEASE CONTACT YOUR LOCAL STERIS REPRESENTATIVE.**

Processing Capacity up to 400 Trays Hours Operation - 24 Peak Processing Hours - 12		
FLOOR PLAN Design Template Sterile Processing Department PROPOSED EQUIPMENT LAYOUT		
 <small>           SYSTEMS PROJECT PLANNING            STERIS Corporation            5900 Healey Road            Mentor, OH 44860 USA            800-548-4837         </small>	Dr: t11 Ck: Date: 23 October 2013 Scale: 3/16" = 1'-0"	H- D- PAGE 1 OF 3

REVISION SCHEDULE		
REVISION #	DESCRIPTION	PLANNER/DATE

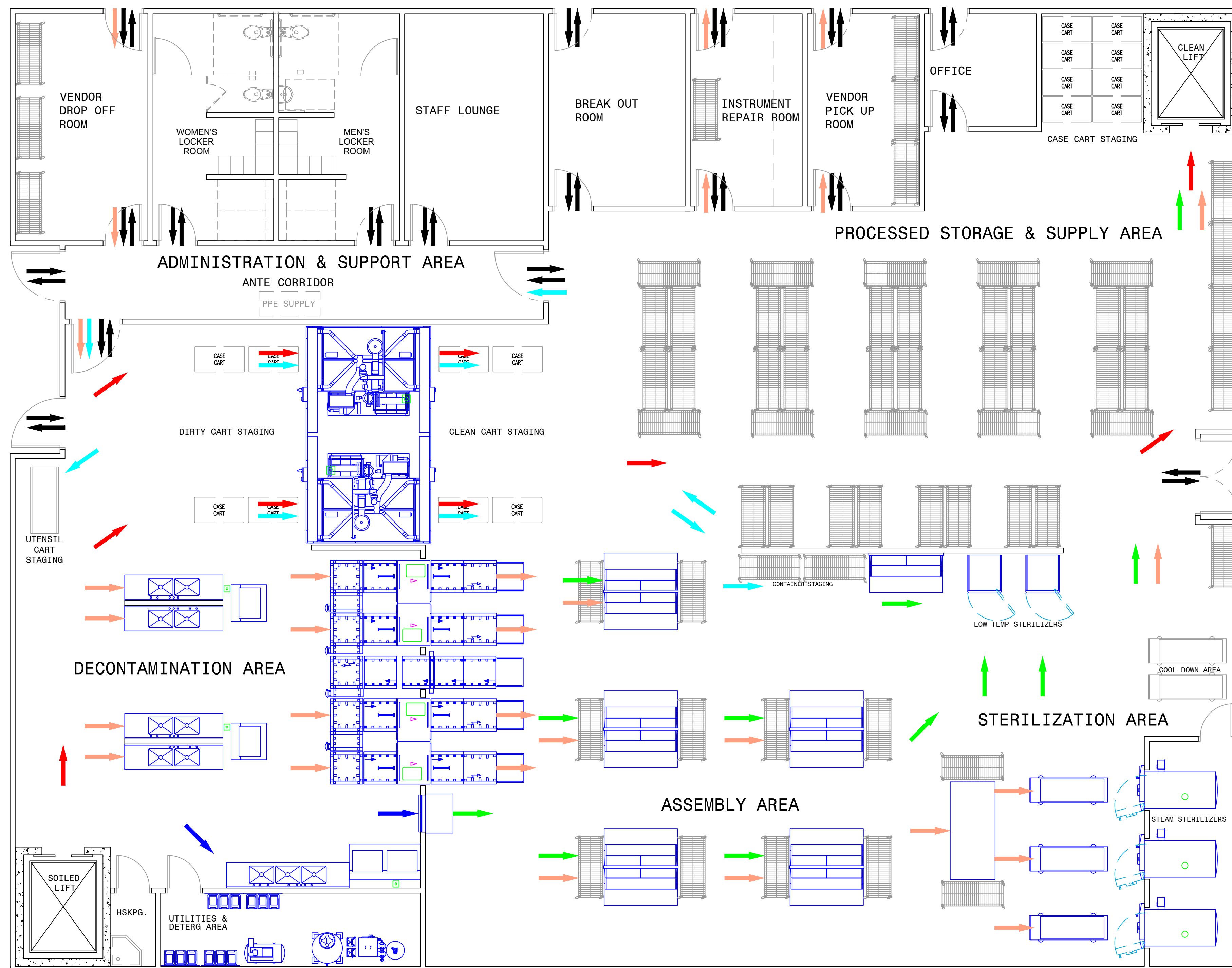


FLOW LEGEND	
CASE CARTS	→
INSTRUMENTS	→
HAND WASH ITEMS	→
LOW TEMP INSTRUMENTS	→
UTENSIL CARTS	→
STAFF	→

## STERILE PROCESSING DEPARTMENTAL FLOW

**\*THIS SAMPLE LAYOUT IS INTENDED TO ENCOURAGE DISCUSSIONS FOR AN INTEGRATED DESIGN PROCESS. FOR A COMPREHENSIVE CAPACITY ASSESSMENT AND OPTIMIZED DESIGN FOR YOUR PROCESSING DEPARTMENT, PLEASE CONTACT YOUR LOCAL STERIS REPRESENTATIVE.**

Processing Capacity up to 400 Trays Hours Operation - 24 Peak Processing Hours - 12		
FLOOR PLAN Design Template Sterile Processing Department PROPOSED EQUIPMENT LAYOUT		
<p>STERIS® SYSTEMS PROJECT PLANNING STERIS Corporation 15900 Healey Road Meriden, CT 06450 USA 800-548-4837</p>	Dr: t11	H-
	Ck:	D-
Date: 23 October 2013	Scale: 3/16" = 1'-0"	PAGE 1 OF 3

REVISION SCHEDULE		
REVISION #	DESCRIPTION	PLANNER/DATE

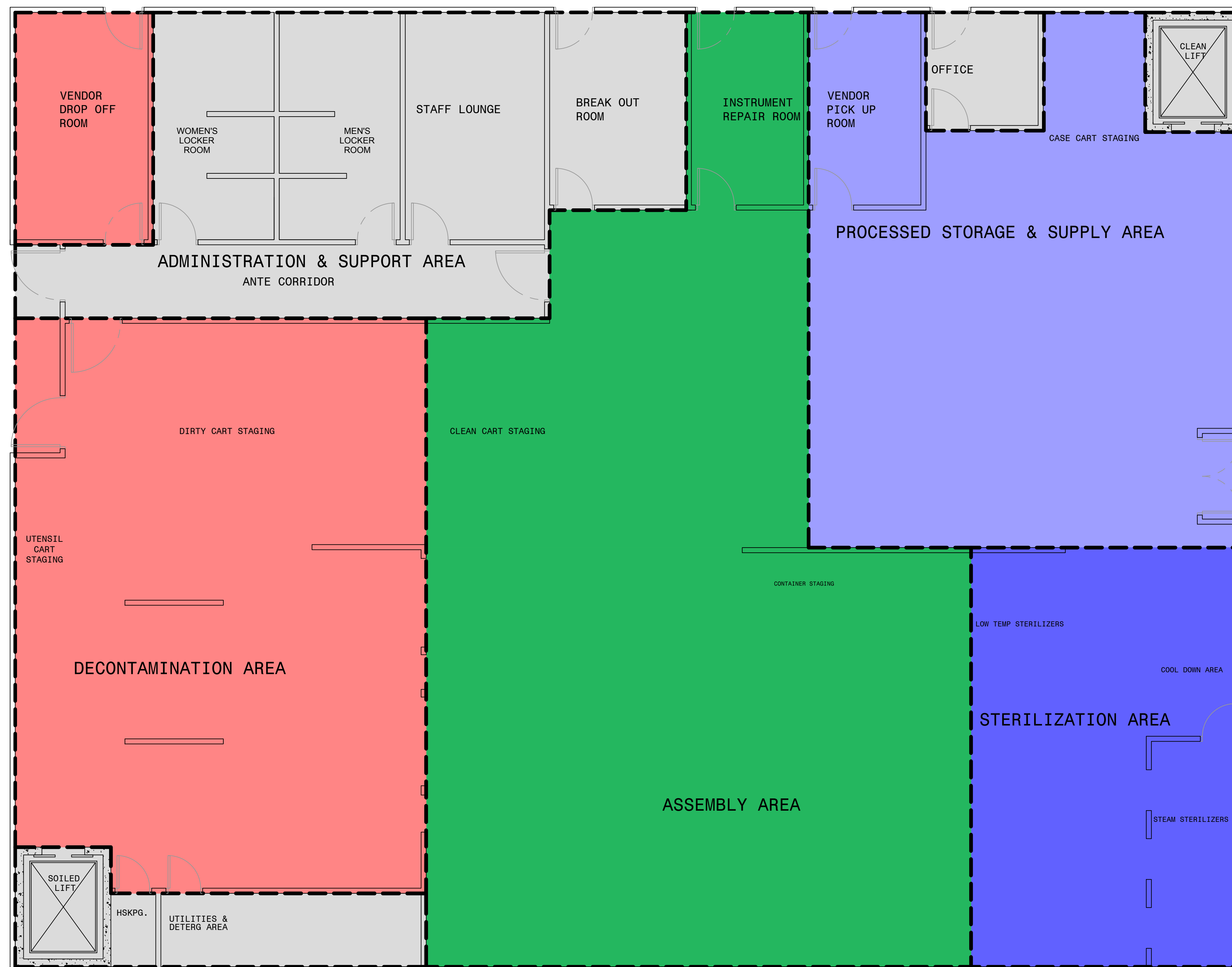

STERILE PROCESSING OVERALL AREA: 7450-8000 ft <sup>2</sup>		(100%)
<span style="color: red;">■</span>	DECONTAMINATION AREA: 1750-1850 ft <sup>2</sup>	(~20%)
<span style="color: green;">■</span>	ASSEMBLY AREA: 2250-2400 ft <sup>2</sup>	(~30%)
<span style="color: blue;">■</span>	STERILIZATION AREA: 750-850 ft <sup>2</sup>	(~10%)
<span style="color: lightblue;">■</span>	PROCESSED STORAGE & SUPPLY AREA: 1300-1400 ft <sup>2</sup>	(~20%)
<span style="color: grey;">■</span>	ADMINISTRATION & SUPPORT AREA: 1400-1500 ft <sup>2</sup>	(~20%)

## STERILE PROCESSING DEPARTMENT AREA DISTRIBUTION

**\*THIS SAMPLE LAYOUT IS INTENDED TO ENCOURAGE DISCUSSIONS FOR AN INTEGRATED DESIGN PROCESS. FOR A COMPREHENSIVE CAPACITY ASSESSMENT AND OPTIMIZED DESIGN FOR YOUR PROCESSING DEPARTMENT, PLEASE CONTACT YOUR LOCAL STERIS REPRESENTATIVE.**

STERIS drawings do not supersede the drawings from the architect of record.

Processing Capacity up to 400 Trays Hours Operation - 24 Peak Processing Hours - 12		
FLOOR PLAN Design Template Sterile Processing Department PROPOSED EQUIPMENT LAYOUT		
 <small>STERIS CORPORATION          5900 Healey Road          Mentor, OH 44860 USA          800-548-4837</small>	Dr: t11	H-
	Date: 23 October 2013	D-
Scale: 3/16" = 1'-0"	PAGE	1 OF 3