


REVISION #	DESCRIPTION	PLANNER/DATE

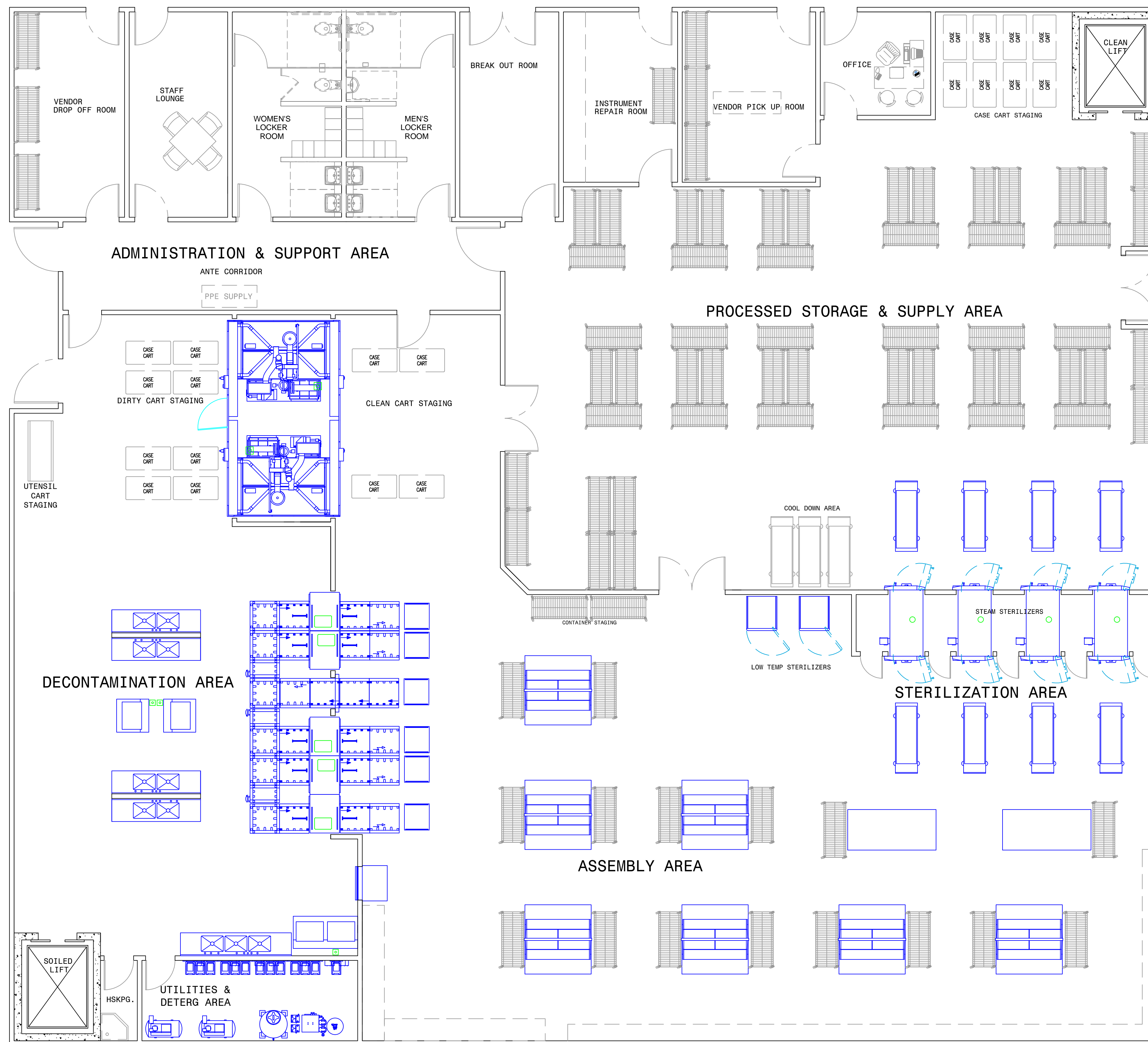

STERILE PROCESSING DEPARTMENT FLOOR PLAN

Processing Capacity up to 600 Trays
 Hours Operation - 24
 Peak Processing Hours - 12


FLOOR PLAN
 Design Template Sterile Processing Department
 PROPOSED EQUIPMENT LAYOUT

***THIS SAMPLE LAYOUT IS INTENDED TO ENCOURAGE DISCUSSIONS FOR AN INTEGRATED DESIGN PROCESS. FOR A COMPREHENSIVE CAPACITY ASSESSMENT AND OPTIMIZED DESIGN FOR YOUR PROCESSING DEPARTMENT, PLEASE CONTACT YOUR LOCAL STERIS REPRESENTATIVE.**

STERIS drawings do not supersede the drawings from the architect of record.

<p>STERIS CORPORATION 3940 HERSHEY ROAD MENTOR, OH 44860 USA 800-548-4837</p>	Dr. t11	H-
	Ck.	D-
	Date 23 October 2013	
	Scale 3/16" = 1'-0"	PAGE 1 OF 3

REVISION #	DESCRIPTION	PLANNER/DATE

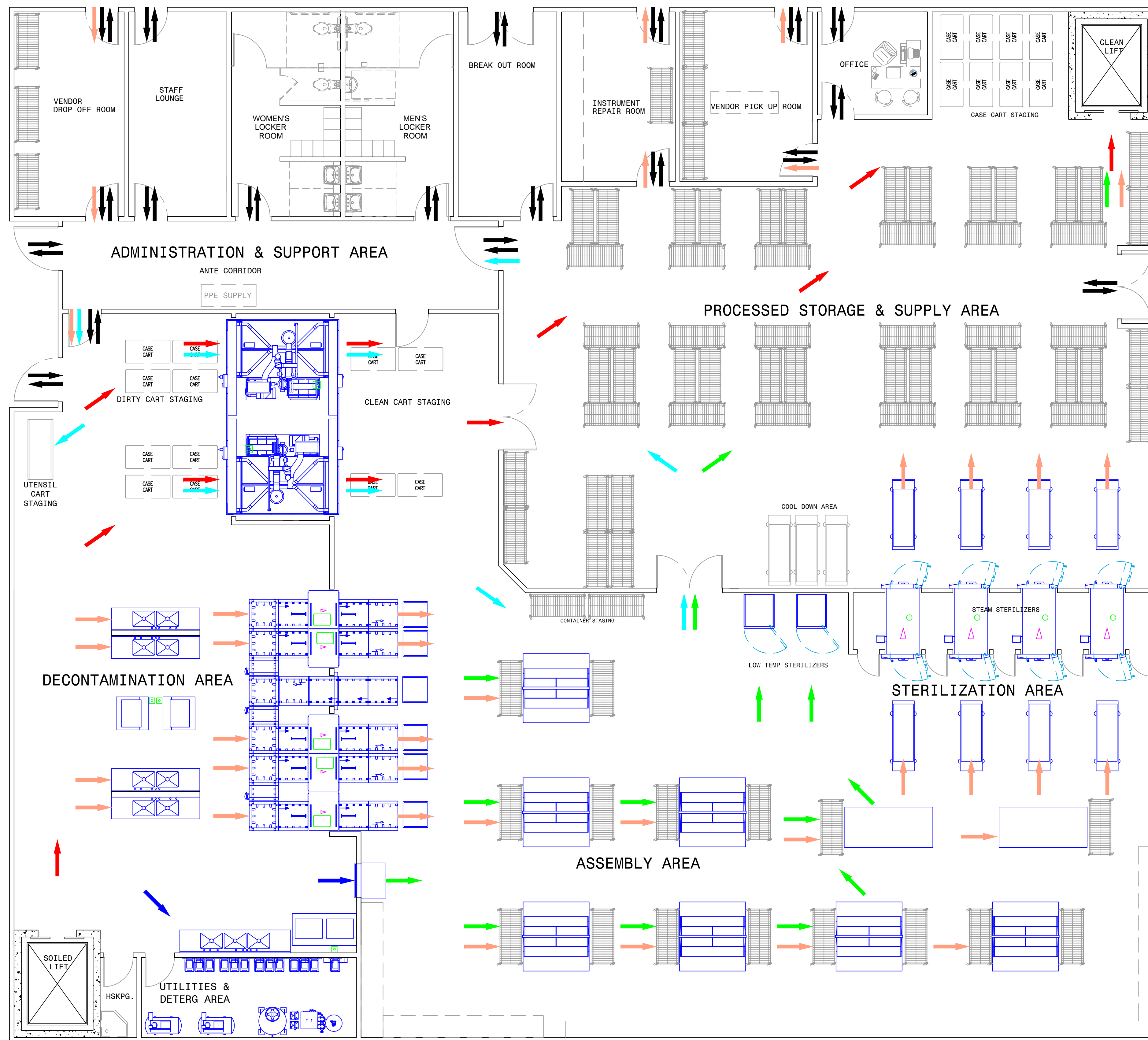

FLOW LEGEND	
CASE CART	→
INSTRUMENTS	→
HAND WASH ITEMS	→
LOW TEMP INSTRUMENTS	→
UTENSIL CARTS	→
STAFF	→


STERILE PROCESSING DEPARTMENTAL FLOW

Processing Capacity up to 600 Trays
 Hours Operation - 24
 Peak Processing Hours - 12

FLOOR PLAN
 Design Template Sterile Processing Department
 PROPOSED EQUIPMENT LAYOUT

***THIS SAMPLE LAYOUT IS INTENDED TO ENCOURAGE DISCUSSIONS FOR AN INTEGRATED DESIGN PROCESS. FOR A COMPREHENSIVE CAPACITY ASSESSMENT AND OPTIMIZED DESIGN FOR YOUR PROCESSING DEPARTMENT, PLEASE CONTACT YOUR LOCAL STERIS REPRESENTATIVE.**

SYSTEMS PROJECT PLANNING STERIS Corporation 5348 Hestley Road Mentor, OH 44860 USA 800-548-4837	Dr. t11	H-
	Ck.	D-
	Date 23 October 2013	PAGE 1 OF 3
	Scale 3/16" = 1'-0"	

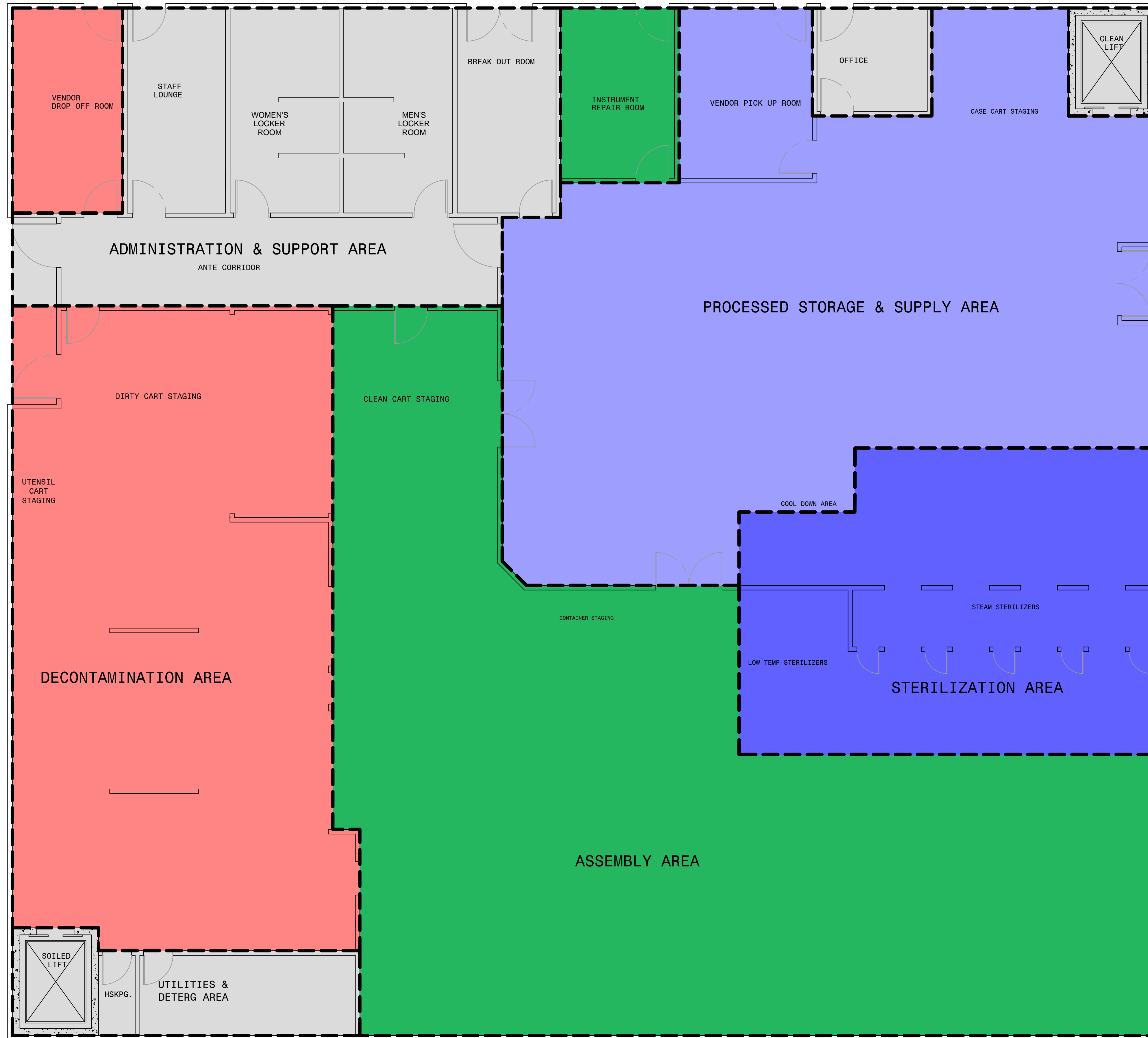

REVISION SCHEDULE		
REVISION #	DESCRIPTION	PLANNER/DATE

	STERILE PROCESSING OVERALL AREA: 9300-9800 ft ²	(100%)
	DECONTAMINATION AREA: 1900-2000 ft ²	(~20%)
	ASSEMBLY AREA: 2700-2800 ft ²	(~30%)
	STERILIZATION AREA: 900-1000 ft ²	(~10%)
	PROCESSED STORAGE & SUPPLY AREA: 1900-2000 ft ²	(~20%)
	ADMINISTRATION & SUPPORT AREA: 1900-2000 ft ²	(~20%)

STERILE PROCESSING DEPARTMENT AREA DISTRIBUTION

***THIS SAMPLE LAYOUT IS INTENDED TO ENCOURAGE DISCUSSIONS FOR AN INTEGRATED DESIGN PROCESS. FOR A COMPREHENSIVE CAPACITY ASSESSMENT AND OPTIMIZED DESIGN FOR YOUR PROCESSING DEPARTMENT, PLEASE CONTACT YOUR LOCAL STERIS REPRESENTATIVE.**

Processing Capacity up to 600 Trays Hours Operation - 24 Peak Processing Hours - 12		
FLOOR PLAN Design Template Sterile Processing Department PROPOSED EQUIPMENT LAYOUT		
 <small>STERIS CORPORATION 3940 HERSHEY ROAD MENTOR, OH 44860 USA 800-548-4837</small>	Dr. t11	H-
	Ck.	D-
	Date 23 October 2013	PAGE
	Scale 3/16" = 1'-0"	1 OF 3